

Drug Profiling Service



Overview


In order to improve anti-cancer drug development and personalized therapies, novel and highly predictive in vitro drug testing systems are needed to examine the characteristics of each individual tumor. Indivumed has developed a preclinical drug testing platform based on precision cut cancer tissue slices

(PCCTS) that is suitable for the investigation of anti-cancer drug effects in the natural tumor microenvironment.

Indivumed offers the cultivation of tumor tissue slices which enables comprehensive drug testing taking into account the complexity of the tumor and its microenvironment.

Workflow

Tumor tissue is collected immediately after resection and subsequently sliced into 500 µm slices. These are randomly transferred to 24-well plates and can be cultivated for up to 48 hours. Due to heterogeneity of the tissue we perform triplicates.

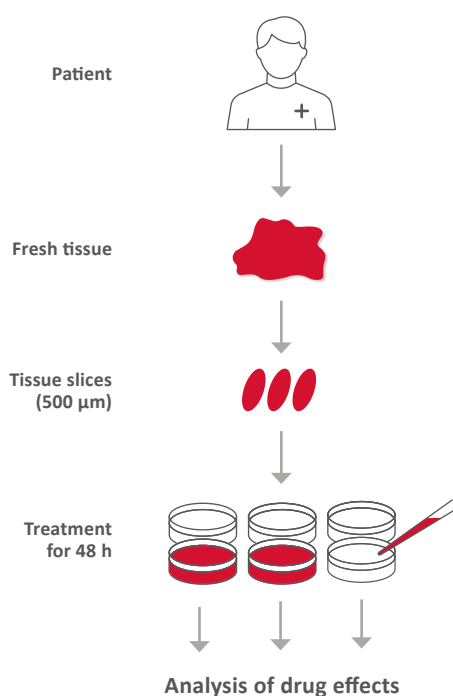

Collection of human tumor tissue material and clinical data

According to customer's specification with regard to tumor entity, pre-treatment, subject group, etc.

Tissue preparation and treatment with compounds

Use Indivumed's tissue expertise for your experiments. Experimental layouts with different compound doses and incubation times up to 48 h.

Analysis of drug effects, for example


- Immunohistochemistry
 - Determination of Tumor Content
 - Immune Infiltration Pattern
 - Co-expression of Biomarkers on a single section
- Gene Expression
 - Mutation Profiling
 - Gene and/or mRNA expression
 - Spatial Transcriptomics
- Multiplex Cytokine Profiling (Meso Scale Discovery)
- Evaluation of Cytotoxicity (LDH-release)

Scientific reporting

Reports include raw data, analysis results, pathology report and clinical information from tissue donors.

Variety of read outs

Evaluation of PD-L1 expression via chromogenic IHC

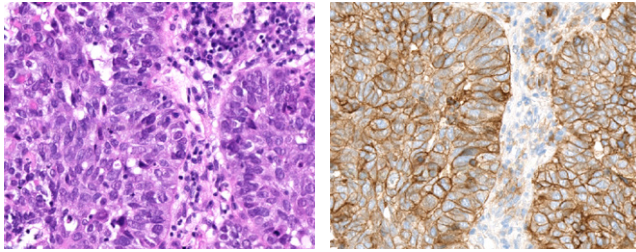


H&E

PD-L1

Anti-PD-L1 IHC of PCCTS at timepoint T0.

PCCTS from a patient diagnosed with NSCLC were analyzed for expression of PD-L1 at T0. Strong membranous anti-PD-L1 staining was detected in the tumor cells.

Identification of T-cell populations via multiplex IHC

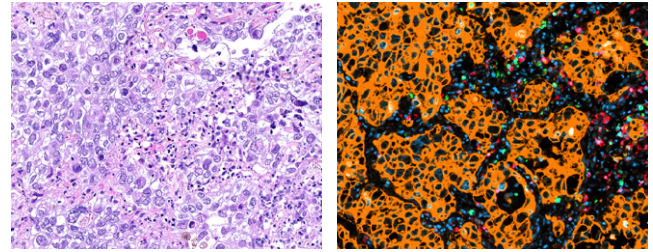


H&E

CD3, CD8, pan-CK

Fluorescent anti-CD8/CD3/pCK mIHC of PCCTS after cultivation for 24 hours.

PCCTS from a patient diagnosed with NSCLC were cultivated for 24 hours and analyzed for expression of CD8 (green), CD3 (red), Pan-CK (orange), DAPI (blue) by mIHC.

Measurement of Cytokine release after Nivolumab treatment

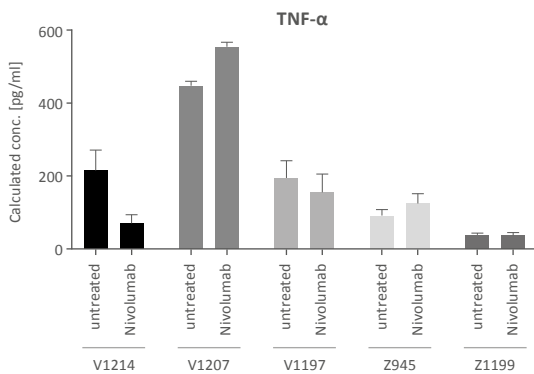

PCCTS from patients diagnosed with NSCLC were treated with 20 µg/ml Nivolumab for 18 hours. Cytokine secretion (TNF-α) was analyzed in supernatants of tissue slices using MSD. Each column represents the mean value with SD in pg/ml of a measurement in duplicates.

Benefits of Indivumed's Precision Cut Cancer Tissue Slice Platform

- Morphological integrity of PCCTS (tumor cells, stromal cells, immune cells) incl. tumor micro-environment
- Complete set of clinical data available from each patient
- Fresh tissue will be delivered within two hours after surgery
- Viable PBMCs can be isolated from the same patient

About Indivumed

Indivumed GmbH is a physician-led, integrated global oncology company for personalized medicine with the world's premier high-content tumor database and highest quality biobank. Indivumed's standard operating procedures are trusted as the global benchmark for biospecimen and clinical data collection.

Our long-standing biobanking expertise is ISO 9001-certified and includes strict protocols for tissue sampling and processing with regard to crucial parameters such as tumor region, ischemia time, and information management.

Contact

Europe Tel.: +49 (40) 41 33 83 18

Email: CustServ@indivumed.com

USA Tel.: +1 (301) 228 9739

www.indivumed.com